

**monogram**<sup>®</sup>  
by beurer

# Instructions



These instructions  
are for the  
Heated  
overblankets



Approved to European Safety Standard EN 60335-2-17












## **Important instructions**

**Please read these instructions carefully  
before you use your heated overblanket,  
and keep them for future use**

## Contents

1. Items included in the package .....3	4.5 Switching on .....5
1.1 Description .....3	4.6 Setting the temperature .....6
2. Important instructions	4.7 Switching off .....6
Retain for future use .....3	4.8 Automatic switch-off .....6
3. Intended use.....5	5. Cleaning and maintenance.....6
4. Operation.....5	6. Storage.....7
4.1 Safety.....5	7. Disposal.....7
4.2 Electrical safety.....5	8. What if there are problems? .....7
4.3 Fitting the heated overblanket to your bad .....5	9. Technical data .....7
4.4 Connecting your heated overblanket .....5	10. Warranty .....8
	11. Helpline.....8

## Explanation of symbols

	Read the instructions!		Machine-washable on an extra-delicate cycle at 30°C
			Do not bleach
	Do not insert pins!		Do not tumble dry
			Do not iron
	Do not use folded or rucked!		Use tumble dryer on the lowest heat only
	Not to be used by very young children (0 - 3 years).		The textiles used for this device meet the stringent human ecological requirements of Oeko Tex Standard 100, as verified by Hohenstein Research Institute.
	<b>WARNING:</b> Warning of risks of injury or health hazards		



**CAUTION:** Safety information about possible damage to appliance/accessories.




**NOTE:** Important information.

## 1. Items included in the package

- 1 Heated overblanket
- 1 Controller
- 1 Instructions for use

### 1.1 Description

- 1. Mains plug
- 2. Mains cable
- 3. Controller
- 4. Indicator lamp
- 5. Sliding switch for ON/OFF and temperature settings
- 6. Plug-in coupling

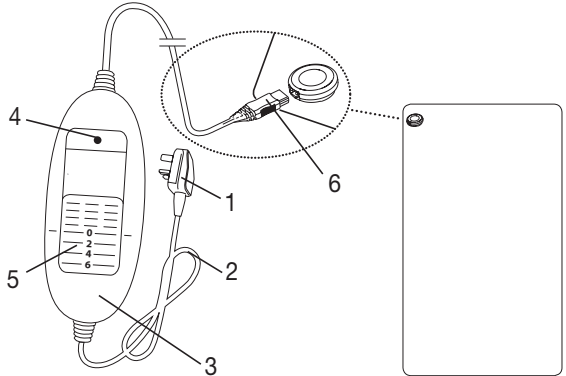

## 2. Important instructions

### Retain for future use



### WARNING

- Non-observance of the following instructions may result in personal injury or material damage (electric shock, skin burns, fire). The following safety and hazard information is not only intended to protect your health and the health of others, it should also protect the product. For this reason, pay attention to these safety notes and include these instructions when handing over the product to others.
- This heated overblanket may not be used by people with reduced sensitivity to heat or by other very vulnerable people as they are unable to respond to overheating.
- This heated overblanket may not be used by very young children (0-3 years old) as they are unable to respond to overheating.
- This heated overblanket may not be used by young children (3-8 years old) unless the controller has been set by a parent or responsible adult and the child has been given adequate instructions on how to use the heated overblanket safely.

- This heated overblanket may be used by children over the age of 8 and by people with reduced physical, sensory or mental skills or a lack of experience or knowledge, provided that they are supervised or have been instructed on how to use the heated overblanket safely, and are fully aware of the consequent risks of use.
- Children must not play with the heated overblanket.
- Cleaning and user maintenance must not be performed by children unless supervised.
- This heated overblanket is not designed for use in hospitals.
- Do not insert needles.
- Do not use when folded or bunched up.
- Do not use if wet.
- This heated overblanket may only be used in conjunction with the controller specified on the label.
- The electrical and magnetic fields emitted by this heated overblanket may interfere with the function of a pacemaker. However, they are still well within the limits: electrical field strength: max. 5000 V/m, magnetic field strength: max. 80 A/m, magnetic flux density: max. 0.1 millitesla. Please consult your doctor and the manufacturer of your pacemaker before using this heated overblanket.
- When this heated overblanket is switched on
  - do not place any objects on it (e.g. suitcase or washing basket)
  - Heat sources, such as hot water bottles, heating pads or similar, must not be placed on the heated overblanket.
- The electronic components in the controller warm up when the heated overblanket is in use.  
For this reason, the controller must never be covered or placed on the heated overblanket when it is in use.
- Do not pull, twist or make sharp bends in the cables.
- Please check this heated overblanket frequently for signs of wear and tear or damage.  
If any such signs are evident, if the heated overblanket has been used incorrectly or if it no longer heats up, it must be checked by the manufacturer before being switched on again.
- If the mains connection cable of the heated overblanket is damaged, it must be replaced by the manufacturer, its customer service department or a similarly qualified person to avoid any risk.
- Always observe the notes on operation (section 4), cleaning and maintenance (section 5) and storage (section 6).
- If you should have any questions about using our devices, please contact Customer Services.

### 3. Intended use



#### CAUTION

This heated overblanket is only designed to warm up the human body. It may only be used as a cover, and not underneath other bedding.

### 4. Operation

#### 4.1 Safety



#### CAUTION

Your heated overblanket has been designed and fully tested to the European Safety Standards. It has been certified by the Intertek Testing and Certification Ltd. UK and shows the safety mark “BEAB Approved”.

This heated overblanket is equipped with a SAFETY SYSTEM. This sensor technology prevents the entire surface of the heated overblanket from overheating by automatically switching off the heated overblanket in the event of a fault.

If the safety system has switched off the heated overblanket, the indicator lamp is not illuminated any more even when the heated overblanket is switched on.

Please note that the heated overblanket is not operable any more due to safety reasons in this case and has to be sent to the service address indicated.

#### 4.2 Electrical safety

Your heated overblanket has a UK plug. If your heated overblanket becomes damaged or worn, you should have the plug top replaced by a qualified electrician who understands how to carry out this repair safely. The plug needs a 3-amp fuse (BS1362 rated).

#### 4.3 Fitting the heated overblanket to your bed

Place the top and bottom sheets on the bed as usual. Spread the heated overblanket onto the bed with the cable entry at the bottom of the bed. The lead(s) should be placed under the bedding and not between the mattress and the base. Place the heated overblanket over the remainder of the bedding.




#### CAUTION

Make sure that the heated overblanket is spread out entirely flat and that it will not rumple up or crease when in use.

#### 4.4 Connecting your heated overblanket

- To operate the heated overblanket connect the controller to the heated overblanket by plugging in the connector.
- Then plug in the supply plug into the mains outlet.

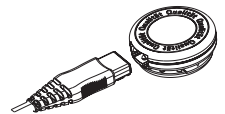

#### 4.5 Switching on

To switch the heated overblanket on, set the sliding switch for ON/OFF and temperature settings to setting 1, 2, 3, 4, 5 or 6.

When it is switched on, the indicator lamp (4) lights up.

## 4.6 Setting the temperature

Level 0: OFF

Level 1: minimum heat

Setting 2 - 5: individual heat setting

Level 6: maximum heat



### WARNING

If the heated overblanket is being used over several hours, we recommend that you set the lowest temperature setting on the controller to avoid overheating the person using the blanket.



### NOTE

The fastest way to warm up the heated overblanket is to initially set the highest temperature setting.

## 4.7 Switching off

To switch off the heated overblanket, slide the sliding switch for ON/OFF and temperature settings to the OFF (0) position.



### NOTE

If the heated overblanket is not being used for several days, slide the sliding switch for ON/OFF and temperature settings to the OFF (0) position and remove the mains plug from the socket. Then disconnect the controller from the heated overblanket by unplugging the plug-in coupling.

## 4.8 Automatic switch-off

This heated overblanket switches off automatically after 12 hours. Once the heated overblanket has switched off, it is no longer heated and the indicator lamp flashes.

To switch on the heated overblanket again, first set the sliding switch for ON/OFF and temperature settings to the OFF (0) position for 5 seconds. After about 5 seconds it is possible to switch it on again.

# 5. Cleaning and maintenance



### WARNING

Before cleaning, always remove the mains plug from the socket first. Then disconnect the controller from the heated overblanket by unplugging the plug-in coupling. Otherwise there is a risk of electric shock.



### CAUTION

- The controller must never come into contact with water or other liquids, as this may cause damage.

To clean the controller, use a dry, lint-free cloth. Do not use any chemical or abrasive cleaning agents. Small marks on the heated overblanket can be removed with a damp cloth or sponge and, if necessary, with a little liquid detergent for delicate laundry.



### **CAUTION**

- Please note that the heated overblanket may not be chemically cleaned, wrung out, tumble dried, put through a mangle or ironed.

Set the washing machine to a particularly gentle wash cycle at 30°C (wool cycle).

Use a delicate laundry detergent and measure it out according to the manufacturer's instructions.



### **CAUTION**

- Please note that frequent washing of the heated overblanket has a negative effect on the product.

The heated overblanket should therefore be washed in a washing machine a maximum of 5 times during its life. Immediately after washing, reshape the heated overblanket to its original dimensions while it is still wet and spread it out flat over a clothes horse to dry.



### **CAUTION**

- Do not use pegs or similar items to attach the heated overblanket to the clothes horse. Failure to do so may damage the heated overblanket.
- Do not reconnect the controller to the heated overblanket until the plug-in connector and heated overblanket are completely dry. Failure to do so may damage the heated overblanket.



### **WARNING**

- Never switch the heated overblanket on to dry it! Otherwise there is a risk of electric shock.

## **6. Storage**

If you do not plan to use the heated overblanket for a long period of time, we recommend that you store it in the original packaging.



### **CAUTION**

Allow the heated overblanket to cool down before folding. Otherwise the heated overblanket may be damaged. Do not crease the heated overblanket by placing items on top of it during storage.

## **7. Disposal**

Please dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment). If you have any queries, please contact the appropriate local authorities.

## 8. What if there are problems?

Problem	Cause	Solution
The indicator lamp is not illuminated while - the controller is properly connected to the heated overblanket - the mains plug is connected to a working socket - the sliding switch for ON/OFF and temperature settings is set to 1, 2, 3, 4, 5 or 6	The safety system has switched the heated overblanket off permanently.	Send the heated overblanket and the controller for servicing.

## 9. Technical data

	Type	Supply Voltage	Input
<b>Overblanket Single</b>	B 10	230–240V / 50–60 Hz	100–110 W
<b>Overblanket Double</b>	B 10	230–240V / 50–60 Hz	2x 100–110 W
<b>Overblanket King Size</b>	B 10	230–240V / 50–60 Hz	2x 100–110 W

## 10. Warranty

Our product is guaranteed for 3 years against material defects or flaws.

This guarantee does not apply to:

- damage resulting from improper use
- parts affected by wear and tear
- faults already known to the customer at the time of purchasing the product
- faults of the customer.

The customer's statutory rights are not affected by the guarantee.

When putting in a claim under guarantee, the customer must provide proof of purchase. Guarantee claims are to be sent to Service-Adress, within a period of 3 years after purchasing the product. Where a claim is justified, the customer has the right to the product being repaired in our own workshop or in a workshop authorised by us.

The customer has no right to further claims (on the basis of the guarantee).

## 11. Helpline

### Distributors:

 Beurer UK Ltd, Suite 10a, Stonecross Place, Yew Tree Way,

 Golborne, WA3 2SH, Tel. 0844 488 9100, Fax 0844 488 9101

### Manufacturer:

Beurer GmbH  
Soeflinger Strasse 218  
89077 Ulm/Germany



Website: [www.beurer.com](http://www.beurer.com)

We have a policy of continuously improving product design and quality.  
We have the right to improve our heated overblanket models at any time.